

## C-Table MKII - Ordering Information

REFERENCE	DESCRIPTION
CTABLE002	C-Table MKII, Imaging Table for Pain Management, Carbon Fiber Top, "Two Motion" 58 - 98 cm and Longitudinal up to 25 cm
CTABLE001*	C-Table MKII, Imaging Table for Pain Management, Carbon Fiber Top, "One Motion" 58 - 98 cm


### OPTIONAL ACCESSORIES

CTABLEAC01	Full Length Accessory Rails (Stainless Steel, Pair)
CTABLEAC02	Lithium Battery Operating System
CTABLEAC03	Handle Bar (Stainless Steel)
CTABLEAC04	Mattress with Face Hole (Visco Foam/Anti-Static)
CTABLEAC05*	Foot Switch Controls (1 Channel)
CTABLEAC06	Foot Switch Controls (2 Channels)

\* An entry-level "One Motion" model remains available. Inquire for the possibilities.

# C-TABLE MKII

## Surprisingly smart & simple!

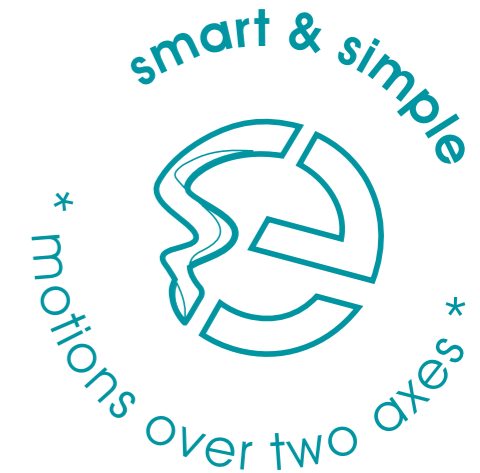

### The Key to Medical Innovation

Innovation in the fields of Anesthesiology and Pain Management is at the heart of Equip Medikey. Joined by true pioneers, we develop new products with unique features, for safer anesthesia and effective pain treatments. Valuable solutions for many practitioners and patients, all over the world.



[www.equip.nl](http://www.equip.nl)

Equip Medikey BV / Edisonstraat 16J / 2809 PB Gouda / The Netherlands

T +31 182 57 32 93 / [info@equip.nl](mailto:info@equip.nl) / [www.equip.nl](http://www.equip.nl)

Rabobank Utrecht / IBAN NL05 RABO 0158 4999 05 / BIC/SWIFT RABONL2U / KvK 24400829 / VAT NL817000240B01



## C-Arm Imaging Table designed for Pain Management


- With a 132 cm unobstructed artefact free imaging area
- Perfectly manoeuvrable and motions over 2 axes
- Lowered to highest height: 58 - 98 cm
- Longitudinal movement of up to 25 cm
- Safe working load of 200 kg
- Wireless controllers

## C-Table MKII for Pain Management

The new C-Table is developed for a wide range of pain management procedures. It combines 132 cm of unobstructed artefact free imaging platform with proven column lift technology, and its low entry height provides easy access for most challenged patients. The line of options and available accessories now also includes wireless hand and foot controllers.

The tables are equipped with a linked brake castor system, operated with a single foot bar. It adjusts easily between three different configurations: fully locked, locked into line (longitudinal movement), and swivel (360 degrees free movement). The head end castors are operated individually. Together, this makes C-Table perfectly manoeuvrable.

Positioning a patient is effortless with this very light treatment table and the additional electric motions over two axes it provides.

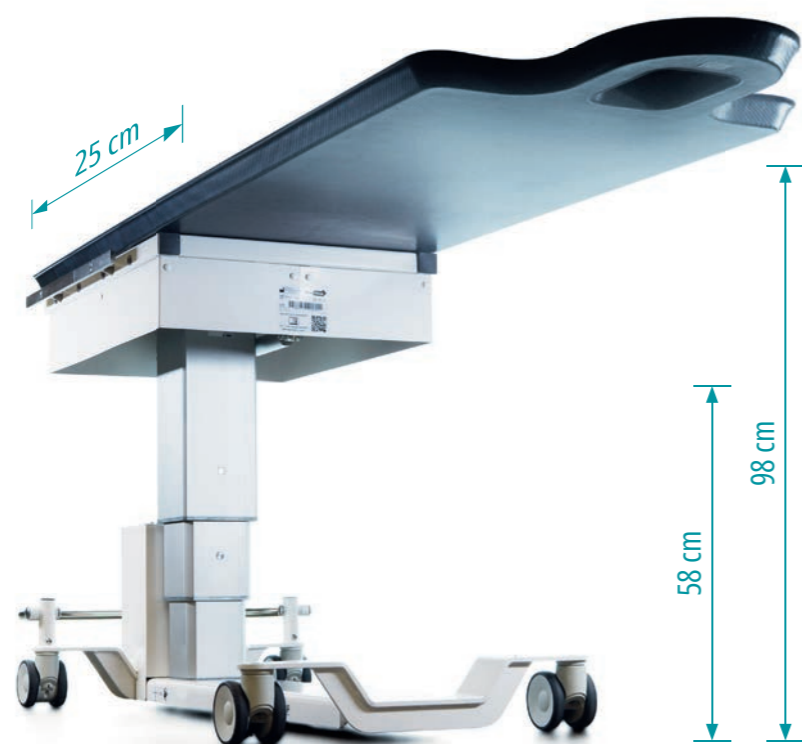

*C-Table MKII offers electric motions over two axes, with a lowered height of 58 cm, a highest height of 98 cm and longitudinal movement of up to 25 cm.*



### C-Table is perfectly manoeuvrable

1. In line (longitudinal movement)
2. Running freely (360° of movement)
3. Fully locked



Safe working load of 200 kg

## STANDARD FEATURES

- Unobstructed artifact free imaging area: 132 cm
- Large carbon fiber top: 208 x 61 cm
- Narrow, lowered base: 12,6 cm
- Lowered height: 58 cm
- Highest height: 98 cm
- Longitudinal movement: up to 25 cm
- Safe working load: 200 kg
- Wireless hand controller
- Rear braking system with easy-to-use foot bar
- Adjustable, medical grade, anti-static castors
- Easy access emergency stop button
- Readily accessible power socket
- Top frame potential grounding point
- CE marked

## AVAILABLE OPTIONS

- Full length accessory rails (pair)
- Stainless steel handle bar
- Battery backup (Lithium)
- Wireless foot controller

## WARRANTY

- Frame: 7 years / Moving parts: 3 years
- Optional battery: 12 months